

Petro Scientific LLC



Density Meter

DM-3005 / DM-3005A / DM-3004

Advanced U-tube vibration

The U-tube vibration densitometer only needs 2~3mL of sample at least, which greatly saves your precious experimental samples.



Peltier temperature control provides you with experimental results at precise temperature.

Excellent display and operation logic

High-definition display of U-tube real-time images, convenient for you to check whether there are bubbles in the U-tube

Follow the built-in SOP to operate, you will get it once you use it.



Support automatic injection- higher throughput

- Standard peristaltic pump to automate your injection, measurement and cleaning processes.
- Ensure accurate results without operator influence.
- Reduce exposure to hazardous materials during injection.




More Features

Complies with ASTM D4052, D5002, D3505, ISO 12185 petroleum standards and all relevant pharmacopoeias.

Four levels of user management permissions, data reporting, electronic signature, and optional audit tracking function.

Built-in network module, equipped with multiple communication interfaces, data can be saved, exported and printed.

<https://www.petro-scientific.com/>

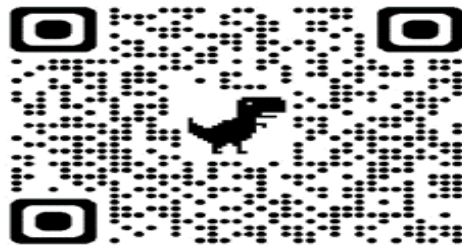

www.petro-scientific.com



Do more with less

Petro Scientific — One-stop Supplier of Scientific Instrument Solutions

Petro Scientific LLC

74-S, Virginia Ave, Penns Grove, NJ 08069 USA

Tel: 856-285-1758

info@petro-scientific.com

www.petro-scientific.com